

Sterican®


Single-use needles for injection  
and aspiration





A complete range of needles with the  
diameters and lengths for every indication


**B | BRAUN**

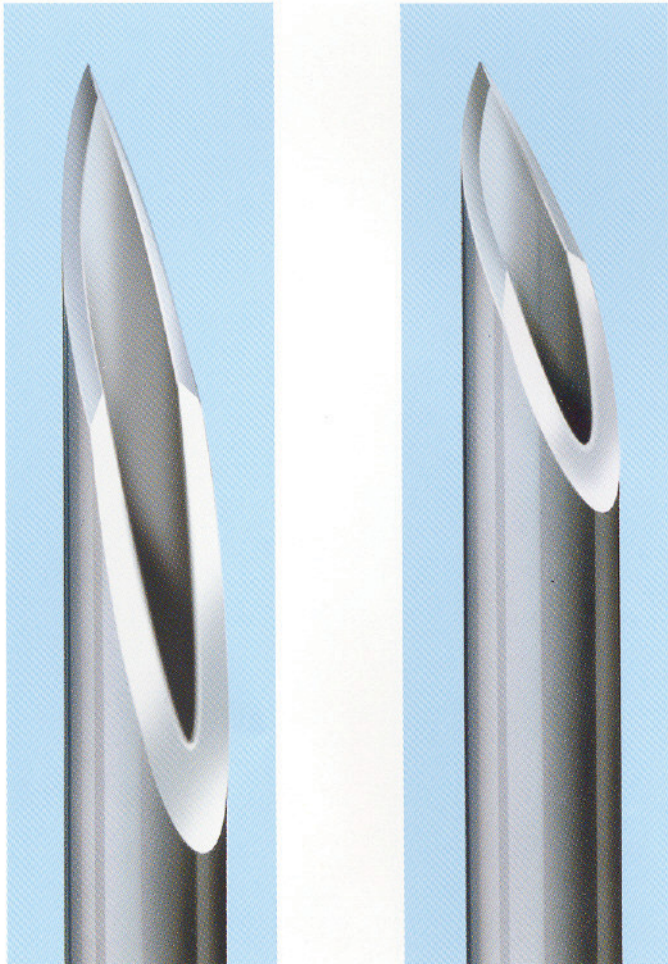Single-use needles for individual indications conforming to ISO 7864/DIN 13097



- 3-facet bevel**
  - for minimal pain upon puncture
- Bevel length according to indication**
  - marked with symbol for easier identification
- Long bevel**  **BL/LB**
  - for easy venipuncture of normal veins
- Short bevel**  **BC/SB**
  - for venipuncture of difficult-to-access veins
- Posterior of bevel eye gently rounded**
  - reduces tissue trauma
- Needle made of stainless chromium nickel steel**
  - with extremely smooth surface and light silicone coating


**Flow rate comparison:**

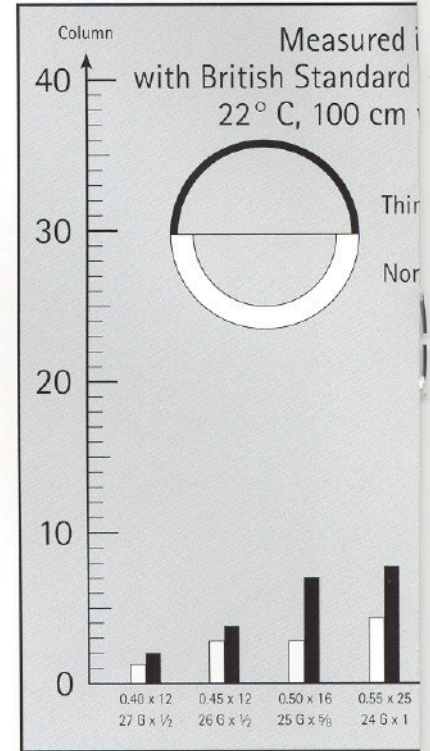

Standard/hypodermic needle range



Special needles for special indications

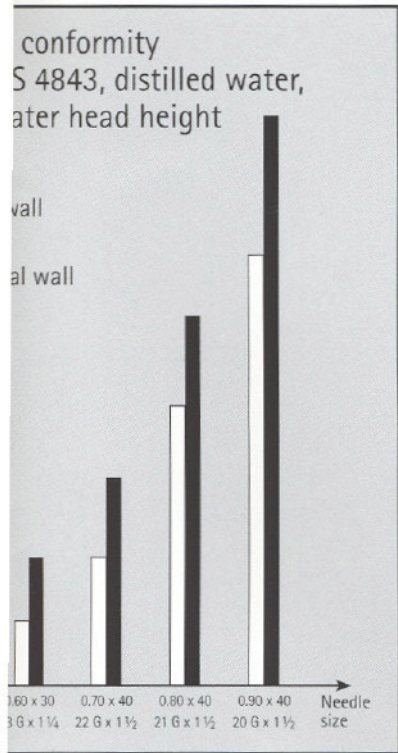

Transparent Luer Lock  
plastic hub made of  
polypropylene

- visualization of blood rapidly shows successful venipuncture

Colour-coded hub

- identifies the outer diameter of the needle, conforms to ISO 6009



Colour-coded dispenser box

- for easier identification and selection of the different sizes
- dispenser boxes containing 100 needles
- bar coded for inventory control



Dental needles

Single-peel packaging

- hygienic removal of sterile needle
- easy-grip peel-off tab



# Sterican®

## Sales packs

Pack size 1000 (10 boxes of 100) Sterican® needles conforming to ISO 7864 and Sterican® Z needles conforming to ISO 7885

Indications	Bevel	Outside diameter mm	Colour code	Length mm	Gauge	Pravaz	Art.-No./Ref.	Pharmacy code		
Sterican® single-use needles Standard range (hypodermic needle) for intramuscular, subcutaneous, intravenous and intraarterial injections	BL/LB	0.90		40	20 G x 1 1/2"	Gr. 1	4657519	2050798		
		0.80		40	21 G x 1 1/2"	Gr. 2	4657527	2050806		
		0.70		30	22 G x 1 1/4"	Gr. 12	4657624	2050812		
		0.60		30	23 G x 1 1/4"	Gr. 14	4657640	2050829		
		0.60		25	23 G x 1"	Gr. 16	4657667	2050835		
		0.55		25	24 G x 1"	Gr. 17	4657675	2050841		
		0.45		25	26 G x 1"	Gr. 18	4657683	2050858		
		0.40		20	27 G x 3/4"	Gr. 20	4657705	2050864		
		for insulin	BL/LB	0.45		12	26 G x 1/2"		4665457	2040676
				0.40		12	27 G x 1/2"		4665406	2499682
0.30				12	30 G x 1/2"		4665300	1410036		
Heparin, tuberculin	BL/LB	0.50		16	25 G x 5/8"		4657853	3200691		
Blood sampling	BL/LB	1.10		40	19 G x 1 1/2"		4657799	3033299		
		1.10		30	19 G x 1 1/4"		4665317	2057949		
		0.90		25	20 G x 1"		4657500	7463246		
		0.80		25	21 G x 1"		4657543	7463186		
		0.70		40	22 G x 1 1/2"		4660021	7463192		
Phlebotomy/Blood sampling Needles with short bevel	BC/SB	1.20		40	18 G x 1 1/2"		4665120	2057955		
		1.10		40	19 G x 1 1/2"		4665112	7516712		
Intramuscular, with short bevel	BC/SB	1.20		50	18 G x 2"		4667123	7463223		
		0.90		50	20 G x 2"		4667093	7463217		
Intramuscular, with long bevel	BL/LB	1.10		50	19 G x 2"		4665511	2057978		
		0.80		50	21 G x 2"		4665503	2057961		
Deep intramuscular	BL/LB	0.90		70	20 G x 2 3/4"		4665791	2057984		
		0.80		80	21 G x 3 1/8"		4665465	0492629		
Neural therapy	BL/LB	0.80		120	21 G x 4 3/4"		4665643	4291448		
		0.60		80	23 G x 3 1/8"		4665635	4291431		
		0.60		60	23 G x 2 3/8"		4665600	2359107		
		0.50		40	25 G x 1 1/2"		9186166	2058009		
Dental upper and lower jaw, conduction anaesthesia infiltration anaesthesia, jaw incisors, upper and lower jaw nerve block anaesthesia	BL/LB	0.50		25	25 G x 1"		9186158	2057990		
		0.40		40	27 G x 1 1/2"		9186182	2058021		
		0.40		25	27 G x 1"		9186174	2058015		
		0.70		50	22 G x 2"		4650018	7516729		
Sterican® special needles	BL/LB	1.20		40	18 G x 1 1/2"		4650034			
		2.10		80	14 G x 3 1/8"		4665473	0492635		
Miscellaneous indications	BC/SB	2.10		80	14 G x 3 1/8"		4665473	0492635		
Sterican® blunt – without bevel for draw up and irrigation	blunt	0.80		22	21 G x 7/8"		9180109	8915986		
		0.40		25	27 G x 1"		9180117	8915992		
Dental – Dental injection needles (cartridge) conforming to ISO 7885	BL/LB	0.30		12	30 G x 1/2"		9324550	4471564		
		0.30		25	30 G x 1"		9324500	2550319		
		0.4		25	27 G x 1"		9324518	2255709		
		0.4		42	27 G x 1 5/8"		9324526	2255690		
		0.5		25	25 G x 1"		9324534	2255684		
		0.5		42	25 G x 1 5/8"		9324542	2255715		

**B | BRAUN**  
HOSPITAL CARE

B. Braun Melsungen AG  
P.O. Box 1120  
D-34209 Melsungen

Tel. +49 (56 61) 71-0  
Fax +49 (56 61) 71-20 44  
www.bbraun.com